


# Movix 4 DReam



see thoroughly,  
ensure safety

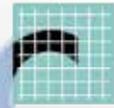

digital mobile x-ray unit

Manoeuvrability  
Efficiency  
Fully digital

**stephani**  
RADIOLOGICAL SOLUTIONS



## O V E R V I E W

The **MOVIX 4 DReam** is a lightweight mobile x-ray unit equipped with a high frequency 4 kW generator and with single or dual **portable Flat Panel Detectors (FPD)** that provides high image quality in any situation. An integrated computer enables instantaneous image display, immediate image processing, local and USB device storage, network archiving and printing.

The lightweight **portable FPD** is as easy to use as cassette. When procedures are completed, the **portable FPD** is placed in the storage compartment and the mobile is ready to move to the next patient.

The **MOVIX 4 DReam** is designed either for outdoor procedures or for in-room examinations.



see thoroughly, ensure safety

## COMPACTNESS



This **MOVIX 4 DReam** is easily foldable for rapid transport in a vehicle or packed into its crate for air shipment. On site, the mobile can be deployed and ready for exposures in as little as **2 minutes**, giving the first image preview within **3 seconds**.

Its manoeuvrability allows all x-ray procedures even in challenging ER and outdoor environments. The **MOVIX 4 DReam** opens new horizons without compromising image quality.

## 3-STEPS READY MOBILE X-RAY UNIT



## E R G O N O M I C S

The dual side hand grips facilitate the tube positioning for efficient access to bedside and polytrauma patients.



The large wheels\* make the **MOVIX 4 DReam** fully dedicated to outdoor situations. This **DR mobile** is designed for operation fields, ER situations and sporting events. With the **MOVIX 4 DReam**, you enter in a new era of digital x-ray.



The wide fork\* increases the stability of the **MOVIX 4 DReam**. Its innovative design facilitates the accessibility where space is tight.



see thoroughly, ensure safety


## ADVANCED TECHNOLOGY

Images are obtained in one single touch with customizable protocols. Preview images are available within **3 seconds** prior to transmission or printout via standardized DICOM interfaces.




The **MOVIX 4 DReam** in association to a **portable FPD** gives you the upper hand in flexibility and versatility, with various features for DICOM data transmission in the latest digital era supports :USB devices, Ethernet cable or wireless.



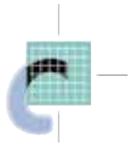

# I M A G E M A N A G E M E N T


The need for an effective approach to information management and storage has never been greater, and STEPHANIX has this approach today. STEPHANIX has surpassed the basic functions of conventional storage to offer unique information management solutions with higher-level features that focus on more efficient management and data preservation.

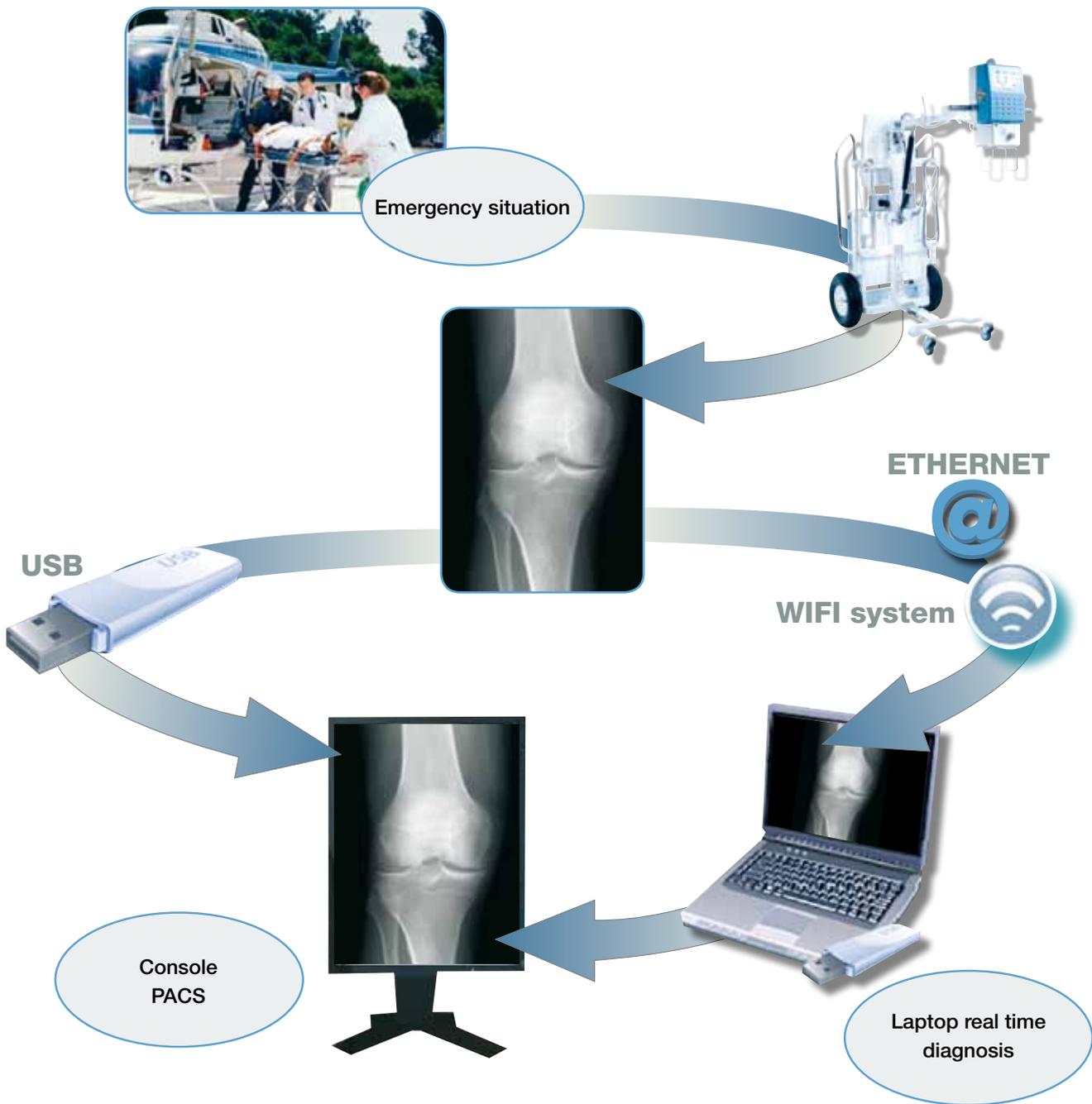

see thoroughly, ensure safety

## A P P L I C A T I O N S

Worldwide shipped with NATO compliance crate, ready in 2 minutes.  
The robust **MOVIX 4 DReam** is made for operations field.

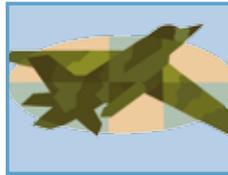

Sport event




The images appear in 3 seconds for real time diagnosis.

The ergonomic tools make you be efficient in any ER situations.

Disaster



NGO

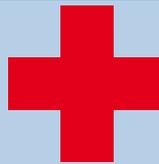

Project in hostile environment has been developed with the Red Cross in French Guyana

As safety comes first, we bring support you to promote **MOVIX 4 DReam** in a special project.

V.I.P



Field event






The **MOVIX DReam** can be easily brought to outdoor events, thanks to the convenience transport in a vehicle.

# Technical specifications

Power	4,0 kW
kV range	from 40 to 115 kV
mA range	from 5 to 100 mA
mAs range	from 0,1 to 200 mAs
Exposure time	from 0,001 s to 10 s
Anatomical protocols	Unlimited
Tube	fixed anode
Focus	0,6 mm and 1,5 mm
Power supply	single phase 230 VCA, 13 A



Detachable solutions

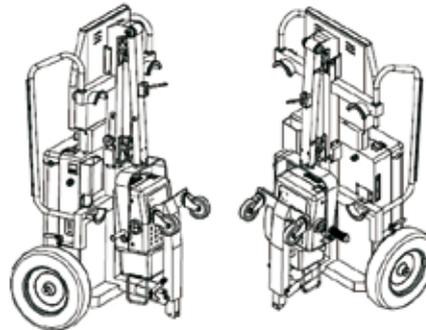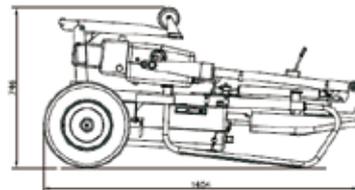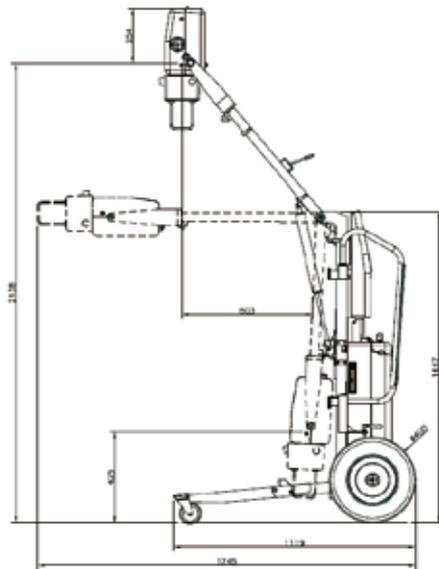

**Head Office**  
 Rue Jean Moulin - ZI du Bayon  
 42150 La Ricamarie - FRANCE  
 Tel : + 33 4 77 47 81 60  
 Fax : + 33 4 77 37 55 19

[www.stephanix.com](http://www.stephanix.com)

**Export department**  
 15, rue du Saule Trapu  
 91882 Massy cedex - FRANCE  
 Tel : + 33 1 60 13 56 25  
 Fax : + 33 4 77 37 55 19